

ARCTIKO LF 650

true cooling specialists



BIOMEDICAL FREEZER

With its low energy consumption, Arctiko's eco-friendly LF series offers flexibility, security, and reliability for users. The Carel controller features an integrated data logger, which can be accessed via Bluetooth.

► INTEGRATED DATA LOGGER

Combining accurate temperature monitoring for optimal sample storage and regulatory compliance, the integrated data logger empowers you with on-the-go control via detailed graphs, alarms, and settings accessible through your phone's Bluetooth connection to the unit's display.

► FORCED AIR COOLING

Forced air cooling provides excellent temperature stability. Forcing the airflow around the stored items ensures quick temperature recovery after door openings.

► KEY LOCK SECURITY


Equipped with a secure key lock, providing an added layer of protection for valuable contents and preventing unauthorized access. This feature ensures the integrity and confidentiality of stored items, making it ideal for environments where security is paramount.

► ENVIRONMENTALLY FRIENDLY REFRIGERANT

This device uses refrigerants with low Ozone Depletion Potential (ODP) and Global Warming Potential (GWP), supporting environmentally conscious operation while maintaining optimal performance and safety standards.



The Carel controller features Bluetooth connection, visual and acoustic alarms, high/low-temperature alarms, open-door alarms and a data logger.

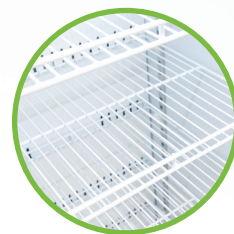

Easily adjustable shelves in order to offer best flexibility for your storage needs.



All models feature durable, cost saving LED lighting for better overview of your samples.



DATA LOGGER



AIR FORCED COOLING



ADJUSTABLE SHELVES






VARIOUS SAFETY FEATURES

ARCTIKO[®] LF 650

true cooling specialists

The ARCTIKO PR, LR and LF series includes the following features, and more:

- Forced Air Cooling
- Low Energy Consumption
- Key Lock on Door
- Carrel controller
- Remote Alarm Interface
- Monoblock Technology
- Eco Friendly Refrigerant

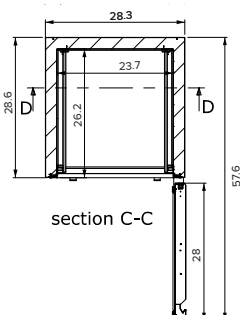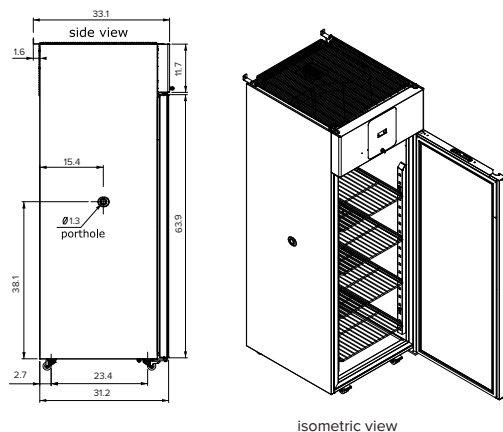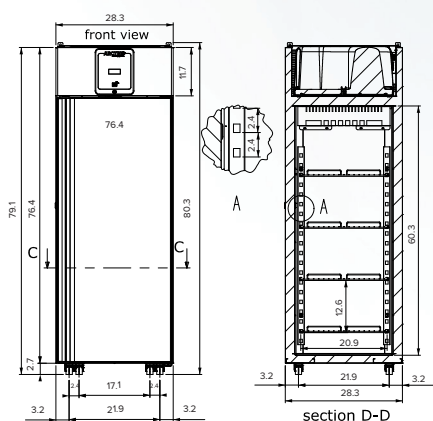

MODEL		LF 650
Main specifications	Item code	1400030-01
	Temperature range (°C)	-25°C to -10°C
	Ambient temperature (°C)	40°C
	Capacity (L)	614
	Capacity (cu.FT)	22
	Exterior dimensions WxDxH (mm)	712x788x2007
	Interior dimensions WxDxH (mm)	610x661x1524
	Weight (KG)	59.87
	Insulation (mm)	60
	Cooling technology	Air Forced
	Controller model / probe type	Carrel iJS / NTC
Power	Power supply (V)	230
	Frequency (Hz)	50
	Power consumption (kWh/24h)	5.43
	AMP (A)	5.2
	Power (Watt)	500
	Horsepower (HP)	1/2
Alarm functions	Alarm (displayed as text, not codes)	Yes
	Visual / acoustic alarm	Yes / Yes
	Power failure alarm	Yes
	Adjustable high / low-temperature alarm	Yes / Yes
	Open door alarm	Yes
	Probe failure alarm	Yes
Data logging & external connection	Low battery alarm	Yes
	Voltage-free contact for remote alarm (e.g. GSM alarm module)	Yes
	Bluetooth data download	Yes
Display features	Electronic data logger	Yes
	Lockable display to prevent accidental shutdowns	Yes
	Temperature Graph Available via app	Bluetooth
Other features	Number of optional reference probes	2
	Automatic defrost	Yes
	Vacuum release port	Yes
	Access port for external probe etc.	Yes
	Light	Yes
	Lock with key	Yes
	Castors (total pcs) / castors with brake (pcs)	4 / 2
	Reversible door(s)	Yes
	Fixed shelves (pcs) / adjustable shelves (pcs)	3
	Max. no. of shelves (pcs) / recommended max. no. of shelves (pcs)	6
	Max load per shelf (lb)	99
	Adjustable drawers (pcs)	Optional
Refrigeration & cabinet	Max. no. of adjustable drawers (pcs)	6
	Max. load per adjustable drawer (lb)	99
	Number of compressors	1
	Refrigerant	R290
	Interior cabinet material	PVC film coated and Anticorrosion treated hot-dip Galvanized steel
	Exterior cabinet material	PVC film coated and Anticorrosion treated hot-dip Galvanized steel
Shipping	Colour (exterior cabinet)	White
	Foam type	Cyclopentane polyurethane
	Shipping weight (kg)	180
Shipping	Shipping dimensions WxDxH (mm)	813x915x2210
	Shipping volume (m³)	1.61