

# ARCTIKO ULUF P510 MV

true cooling specialists



**ULTRA LOW TEMPERATURE | HIGHLY RELIABLE ULT FREEZER WITH NATURAL REFRIGERANT AND LOW POWER CONSUMPTION**

The freezer comes with the extremely powerful independent **TRUE DUAL™** technology that keeps the samples protected at -70 °C if an unexpected failure occurred keeping you worry-free and the samples secured and safely guarded.



The first-class V700 controller is equipped with all necessary alarms and features on top of everything, we developed ECO mode to reduce energy consumption while maintaining optimum uniformity for the stored materials.



NATURAL REFRIGERANT

**REACHES -80°C IN LESS THAN 4 HOURS**

## ▶ RELIABLE FREEZING

This model features a swift pull-down time of about 3 hours, showcasing a strong cooling system with robust insulation from a mix of polyurethane and vacuum-insulated panels.

Even with frequent door openings, it quickly reestablishes a uniform temperature, maintaining a reliable environment unaffected by external temperatures, thanks to its excellent temperature recovery capabilities.

1 See below graph under "Improved recovery time". Ambient temperature 32°C Freezer with no-load.

## ▶ ENVIRONMENTALLY FRIENDLY

Our natural refrigerants boast ODP1 and GWP2 values under 4, featuring CFC/HCFC-free formulations that protect your samples effectively while safeguarding the Earth.

These eco-friendly refrigerants also reduce energy consumption, enhancing the cooling system's cost-effectiveness and maintaining high performance through strong thermodynamic properties.

1 See below graph under "Improved recovery time". Ambient temperature 32°C Freezer with no-load.

## ▶ ULTIMATE SAMPLE PROTECTION

This model incorporates advanced technologies and the very best components for exceptional performance.

Our **TRUE DUAL™** technology prevents sample loss by keeping the freezer operational in the unlikely event of one compressor failing, maintaining -70°C. This ensures your research or samples are protected around the clock, avoiding the risks associated with inferior laboratory equipment.

## ▶ MULTI-VOLTAGE

Our multi-voltage freezers operate between 100 to 240 volts, seamlessly integrating into global settings without additional adapters or transformers.

This adaptability also ensures resilience against power fluctuations, providing safe and stable storage for sensitive materials.



ULUF P510 MV



We designed the model with easily accessible compartments including four removable inner doors that are individually insulated with magnetic catches. This smart function also guarantees easy defrosting and effortless maintenance of the unit.



Ergonomic door handle enables one-handed door operations. Your security will be doubled thanks to both the integrated lock inside the handle and the RFID smart lock Card Key.



The modernized shelves in the new Integraline are perforated and easily adjustable, offering the best flexibility for your storage needs. Extra durable shelves with up to 50 kg.



TRUE DUAL™ TECHNOLOGY



ONEGRIP HANDLE





DIRECT COOLING




VIP VACUUM INSULATED PANELS

\* The pictures shown are for reference only, actual product size may vary. Please read the dimensions of the product in the specifications.

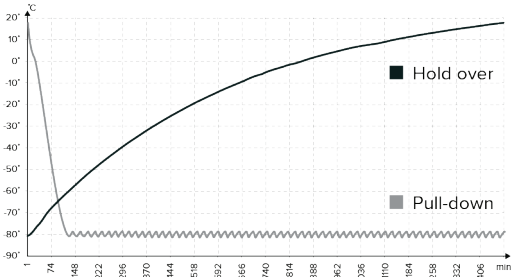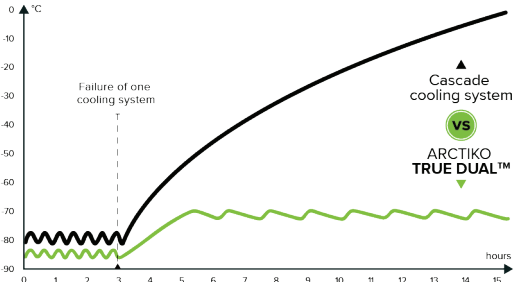

### ▲ IMPROVED RECOVERY TIMES


Pull-down times measures the time taken to reach the correct temperature, the faster the better. Hold over indication the freezers ability to keep the samples safe.

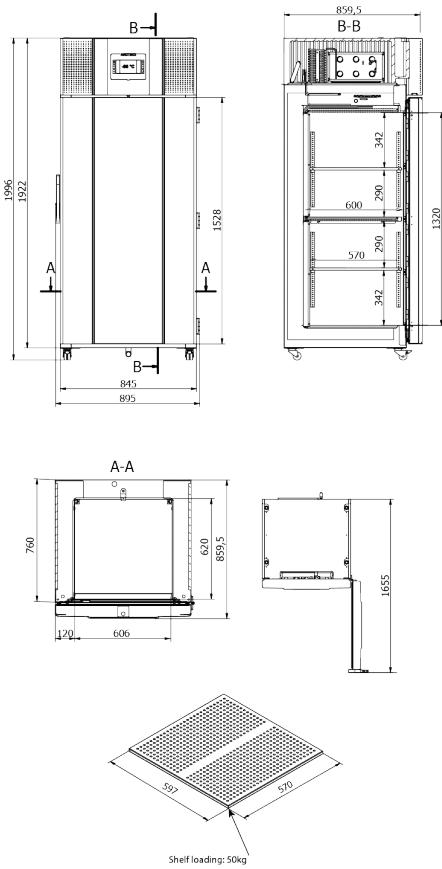

### ◀ CASCADE VS ARCTIKO TRUE DUAL™ System

With a Cascade system (marked black ●) the temperature will rise gradually and the samples will be destroyed.

ARCTIKO's TRUE DUAL™ cooling system (marked green ●) will serve as a backup and the inside temperature will stay at a low temperatures

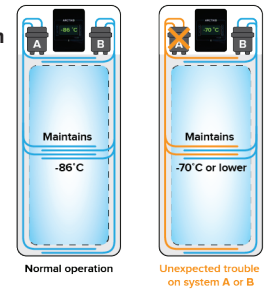

MODEL		ULUF P510 MV
Main specifications	Item code	1000048
	Temperature range (°C)	-40 / -86
	Max. ambient temperature (°C)	32
	Capacity (L)	500
	Exterior dimensions WxDxH (mm)	895x859,5x1999
	Interior dimensions WxDxH (mm)	606x620x1320
	Weight (KG)	233
	Insulation (mm)	120
	Cooling technology	dual system / static
	Controller model / probe type	V700 / PT100
Power	Power supply (V)	100 to 240
	Frequency (Hz)	50/60
	Power consumption (kWh/24h)	8,1
	AMP (A)	-
	Power (Watt)	-
	Horsepower (HP)	-
Alarm functions	Alarm (displayed as text, not codes)	Yes
	Visual / acoustic alarm	Yes / Yes
	Power failure alarm	Yes
	Adjustable high / low-temperature alarm	Yes / Yes
	Open door alarm	Yes
	Probe failure alarm	Yes
	Low battery alarm	Yes
	Battery backup for alarms (approx. hrs.)	72
Data logging & external connection	Voltage-free contact for remote alarm (e.g. GSM alarm module)	Yes
	Electronic data logger	Yes
	RS 485 port	Optional
	USB port for software update & data download	Yes
	Temperature chart recorder	Optional (external mounting)
Display features	Display type	7" touch screen
	Number of optional reference probes	0 (1 reference probe included as standard)
	Password protection for turning unit on / off	yes
	Temperature graph on display	yes
Other features	Vacuum release port	yes
	Antifreeze function	yes
	Access port for external probe etc. (int. diameter mm)	12
	Light	yes
	Lock with key	Yes
	Castors (total pcs) / castors with brake (pcs)	4 / 2
	Sub lids / inner doors (pcs)	4
	Fixed shelves (pcs) / adjustable shelves (pcs)	1 / 2
	Max. no. of shelves (pcs) / recommended max. no. of shelves (pcs)	5 / 5
	Max. load per shelf (kg)	50
	Total capacity of 50 mm boxes 9x9 dividers (pcs)	352
Total capacity of 75 mm boxes 9x9 dividers (pcs)	224	
Total capacity of 96 mm boxes 9x9 dividers (pcs)	160	
Max. capacity of samples in Arctiko 10x10 boxes (pcs)	35,200	
Refrigeration & cabinet	Number of compressors	2
	Refrigerant	GG20
	Interior cabinet material	Stainless steel
	Exterior cabinet material	Painted steel
	Colour (exterior cabinet)	White
	Foam type	Polyurethane
Shipping	VIP vacuum insulated panels	Yes
	Shipping weight (kg)	366
	Shipping dimensions WxDxH (cm)	101x127x225
Shipping volume (m <sup>3</sup> )	2,89	



IEC 61010 / EN IEC 61010 / UL 61010 / CSA C22.2 no 61010