



#### **BIOMEDICAL REFRIGERATOR**

With its low energy consumption, Arctiko's eco-friendly LR series offers flexibility, security, and reliability for users. The Carel controller features an integrated data logger, which can be accessed via Bluetooth.

#### ENERGY STAR CERTIFIED

Energy Star certified, guaranteeing exceptional energy efficiency and reduced environmental impact, while maintaining superior performance in temperaturesensitive storage applications.

## INTEGRATED DATA LOGGER

Combining accurate temperature monitoring for optimal sample storage and regulatory compliance, the integrated data logger empowers you with on-the-go control via detailed graphs, alarms, and settings accessible through your phone's Bluetooth connection to the unit's display.

### FORCED AIR COOLING

Forced air cooling provides excellent temperature stability. Forcing the airflow around the stored items ensures quick temperature recovery after door openings.

# **KEY LOCK SECURITY**

Equipped with a secure key lock, providing an added layer of protection for valuable contents and preventing unauthorized access. This feature ensures the integrity and confidentiality of stored items, making it ideal for environments where security is

### **SNAP/EPA COMPLIANT** REFRIGERANT

Utilizes refrigerants that are fully compliant with the Environmental Protection Agency (EPA) and the Significant New Alternatives Policy (SNAP) regulations, ensuring environmentally responsible operation without compromising performance or safety.



RCTIKE



DATA LOGGER



ADJUSTABLE



SHELVES



VARIOUS



The Carel controller features Bluetooth connection, visual and acoustic alarms. high/low-temperature alarms, open-door alarms and a data logger.



The shelves are easily adjustable to offer the best flexibility for your storage needs




This model features durable cost saving LED lighting for better overview of your samples


\* The pictures shown are for reference only, actual product size may vary. Please read the dimensions of the product in the specifications.





The ARCTIKO PR, LR and LF series includes the following features, and more:


- Forced Air Cooling
- Low Energy Consumption
- Key Lock on Door
- Carel controller
- Remote Alarm Interface
- Monoblock Technology
- Eco Friendly Refrigerant

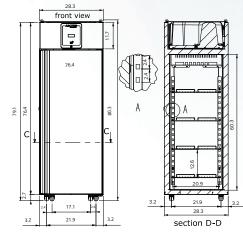

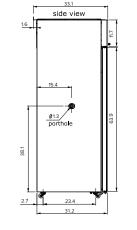

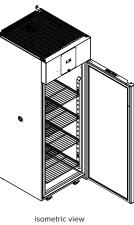

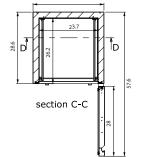

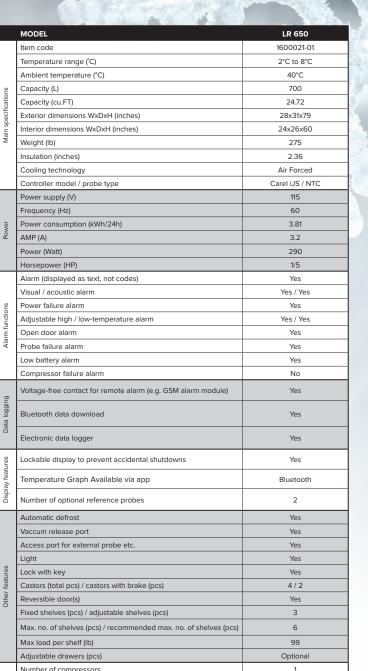

Number of compressors Refrigerant Refrigeration & cabinet Interior cabinet material Anticorrosion treated hot-dip Anticorrosion treated hot-dip Exterior cabinet material Colour (exterior cabinet) Foam type Cyclopentane polyurethane Shipping weight (lb) Shipping dimensions WxDxH (inches) Shipping volume (cu.FT)

ENERGY STA Intertek

R290

PVC film coated and

galvanized steel

PVC film coated and

galvanized steel

White

330

32x36x87

56.86

MAIN DIMENSIONS