

true cooling specialists



BIOMEDICAL REFRIGERATOR

With its low energy consumption, Arctiko's eco-friendly PR series offers flexibility, security, and reliability for users. The Carel controller features an integrated data logger, which can be accessed via Bluetooth.

ENERGY STAR CERTIFIED

Energy Star certified, guaranteeing exceptional energy efficiency and reduced environmental impact, while maintaining superior performance in temperaturesensitive storage applications.

▶ INTEGRATED DATA LOGGER

Combining accurate temperature monitoring for optimal sample storage and regulatory compliance, the integrated data logger empowers you with on-the-go control via detailed graphs, alarms, and settings accessible through your phone's Bluetooth connection to the unit's display.

▶ FORCED AIR COOLING

Forced air cooling provides excellent temperature stability. Forcing the airflow around the stored items ensures quick temperature recovery after door openings.

KEY LOCK SECURITY

Equipped with a secure key lock, providing an added layer of protection for valuable contents and preventing unauthorized access. This feature ensures the integrity and confidentiality of stored items, making it ideal for environments where security is paramount

► SNAP/EPA COMPLIANT REFRIGERANT

Utilizes refrigerants that are fully compliant with the Environmental Protection Agency (EPA) and the Significant New Alternatives Policy (SNAP) regulations, ensuring environmentally responsible operation without compromising performance or safety.





DATA LOGGER



AIR FORCED COOLING



ADJUSTABLE SHELVES





The Carel controller features Bluetooth connection, visual and acoustic alarms, high/low-temperature alarms, open-door alarms and a data logger.



Easily adjustable shelves in order to offer best flexibility for your storage needs.




All models feature durable, cost saving LED lighting for better overview of your samples

^{*} The pictures shown are for reference only, actual product size may vary. Please read the dimensions of the product in the specifications.



The ARCTIKO PR, LR and LF series includes the following features, and more:

- Forced Air Cooling
- Low Energy Consumption
- Key Lock on Door
- Carel controller
- Remote Alarm Interface
- Monoblock Technology
- Eco Friendly Refrigerant

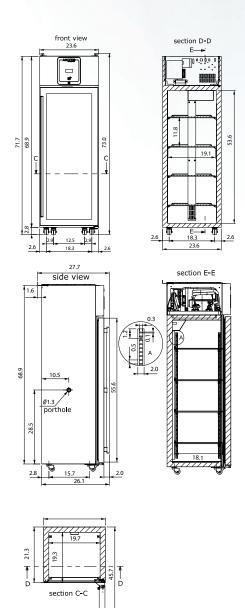

Main specifications	Item code	1700027-01
	Temperature range (°C)	+2 to +8 °C
	Ambient temperature (°C)	40 °C
	Capacity (L)	400
	Capacity (cu.FT)	14.13
	Exterior dimensions WxDxH (inches)	24x24x72
	Interior dimensions WxDxH (inches)	20x19x54
	Weight (lb)	242
	Insulation (inches)	1.97
	Cooling technology	Air Forced
	Controller model / probe type	Carel iJS / NTC
Power	Power supply (V)	115
	Frequency (Hz)	60
	Power consumption (kWh/24h)	3.04
	AMP (A)	3
	Power (Watt)	260
	Horsepower (HP)	1/4
Alarm functions	Visual / acoustic alarm	Yes / Yes
	Power failure alarm	Yes
	Adjustable high / low-temperature alarm	Yes / Yes
	Open door alarm	Yes
	Probe failure alarm	Yes
	Low battery alarm	Yes
Data logging & external connection	Voltage-free contact for remote alarm (e.g. GSM alarm module)	Yes
	Bluetooth data download	Yes
	Electronic data logger	Yes
Display	Lockable display to prevent accidental shutdowns	Yes
	Temperature Graph Available via app	Bluetooth
	Number of optional reference probes	2
Other features	Automatic defrost	Yes
	Access port for external probe etc.	Yes
	Light	Yes
	Lock with key	Yes
	Castors (total pcs) / castors with brake (pcs)	4/2
	Fixed shelves (pcs) / adjustable shelves (pcs)	3
	Max. no. of shelves (pcs) / recommended max. no. of shelves (pcs)	6
	Max Load Per Shelf (lb)	66
Refrigeration & cabinet	Number of compressors	1
	Refrigerant	R290
	Interior cabinet material	Hot-Dipped Galvanized Steel, Anticorrosion Treated, PVC Film Coated
	Exterior cabinet material	Hot-Dipped Galvanized Steel, Anticorrosion Treated, PVC Film Coated
	Colour (exterior cabinet)	White
	Foam type	Cyclopentane polyurethane
Shipping	Shipping weight (lb)	286
	Shipping dimensions WxDxH (inches)	27x28x80
	Shipping volume (cu ET)	34.26

Shipping volume (cu.FT)





34.26